


CEKA ATTACHMENTS PRECI-LINE


CEKA 334



**PLENTY OF SPACE
FOR AN ATTACHMENT?**

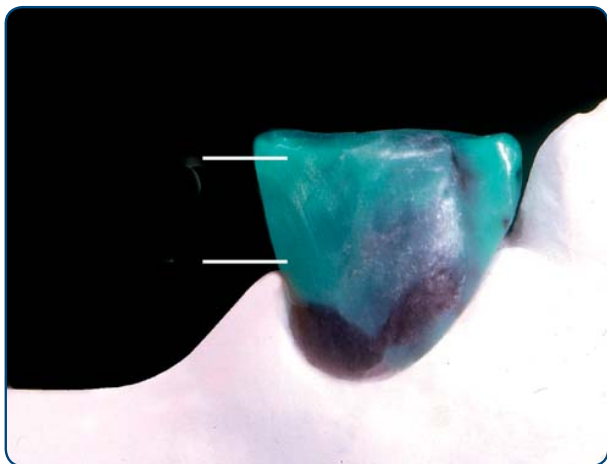

**OF COURSE,
WITH CEKA REVAX!**

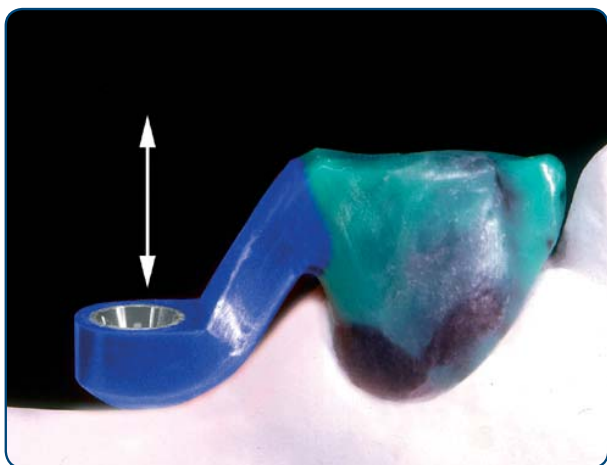

FAMILIAR COMPONENTS

THE PAST



OL 334 NX

2 attachments/package

THE FUTURE



OL 0234 NPR

1 attachment/package

ESTHETICS, STABILITY AND ORAL HYGIENE

TRY IT!



Press the cast-on precision female with the colour-coded side down into the keeper.



Reduce the plastic keeper at the buccal. This creates a natural transition to the crown and more space for esthetics.



Position the profile on the crest with passive contact. Cast in a non-precious alloy.



The attachment female provides adequate space for good contour and esthetics for a natural appearing removable prosthesis.



SWISS MADE



CEKA

CEKA SOL
PRECI-HORIX
PRECI-VERTIX
PRECI-VERTIX AT
PRECI-BAR
PRECI-BALL
PRECI-CLIP
PRECI-SAGIX
PRECI-CLIX AXIAL
PRECI-CLIX EXTRACORONAL
PRECI-CLIX RADICULAR
PRECI-CLIX RADICULAR RC
PRECI-TUBIX/MORTIX
PRECI-POST
PRECI-PROFILE
PERMA-RET
3C-BOND
EXPANDO
CEKA-MAX

CEKA 334