

CEKA ATTACHMENTS PRECI-LINE

PRECI-HORIX®



PRECI-HORIX

SIMPLE AND ECONOMICAL

The PRECI-HORIX bar and rider system, developed by Helmut Hader, is the most popular of all bar systems due to its economy and simplicity.

The excellent retention of the plastic clips, the patient comfort, the ease of processing and durability make PRECI-HORIX a first choice by experts for over 30 years.




- Easily replaceable female riders in three retention levels to choose from
- The bar pattern may be cast in alloy of choice.
- Extremely suitable for use on implants
- Direct support and rotational movement
- The soft snap of the female riders offers the patient security.



PRECI-HORIX HOUSING

SIMPLIFIES CLIP PLACEMENT

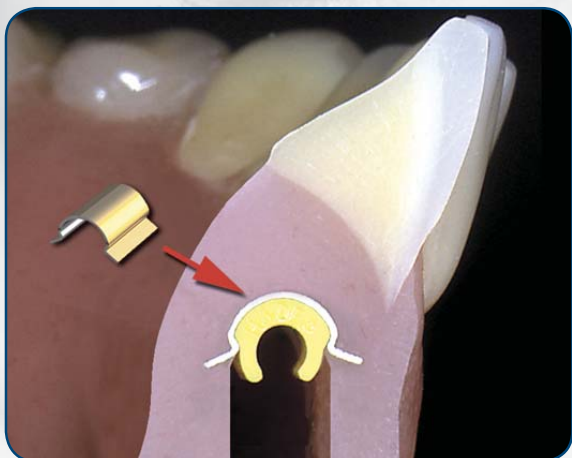




Recommended use is processing the metal housings into the denture simplifying clip placement and then using one of the three retention grips of plastic PRECI-HORIX clips in the metal housings.


- White: reduced retention
- Yellow: normal retention
- Red: increased retention

The prevailing 1.8 mm bar pattern is compatible with PRECI-CLIP riders (see separate leaflet).



EASY PROCESSING

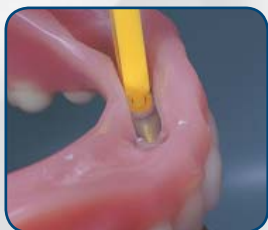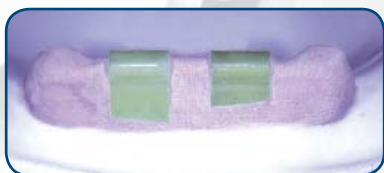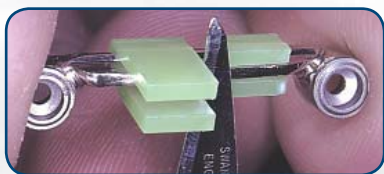

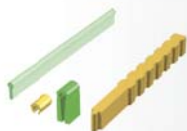


PRECI-HORIX STANDARD

1701

KIT

- 2 plastic bar patterns
- 6 yellow riders
- 6 processing riders
- 1 insertion tool

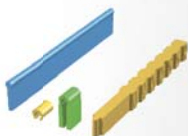

PRECI-HORIX PLASTICWAX

1731

KIT

- 2 PLASTICWAX bar patterns
- 6 yellow riders
- 6 processing riders
- 1 insertion tool



PRECI-HORIX COMBI

1741

KIT

- 1 plastic bar pattern
- 2 yellow riders
- 2 processing riders
- 2 stainless steel housings
- 1 insertion tool



PRECI-HORIX RIDER

normal retention:

1802 yellow 6 pcs

1802B yellow 50 pcs

reduced retention:

1805 white 6 pcs

1805B white 50 pcs

increased retention:

1806 red 6 pcs

1806B red 50 pcs

combination package:

1807 2 + 2 + 2 pcs



PRECI-HORIX METAL HOUSING

Ultrathin stainless steel housing with titanium/nitrate layer (gold-coloured) for a perfect fit of the riders in the prosthesis

1706 6 pcs

1706B 50 pcs



PRECI-HORIX PLASTIC PATTERN

Burn-out plastic pattern.

Cast in any hard alloy.

Length: 50 mm – Height: 4 mm – \varnothing 1.8 mm

1704 6 pcs

1704B 50 pcs



PRECI-HORIX PLASTICWAX PATTERN

Ideal for use with fast set investments.

Cast in any hard alloy.

Length: 50 mm – Height: 8.25 mm – \varnothing 1.8 mm

1734 6 pcs

1734B 50 pcs

WWW.CEKA-PRECI-LINE.COM

SWISS MADE



CEKA REVAX

CEKA SOL

PRECI-HORIX

PRECI-VERTIX

PRECI-VERTIX AT

PRECI-BAR

PRECI-BALL

PRECI-CLIP

PRECI-SAGIX

PRECI-CLIX AXIAL

PRECI-CLIX EXTRACORONAL

PRECI-CLIX RADICULAR

PRECI-CLIX RADICULAR RC

PRECI-TUBIX/MORTIX

PRECI-POST

PRECI-PROFILE

PERMA-RET

3C-BOND

EXPANDO

CEKA-MAX

PRECI-HORIX®

Prospectus no 273.E-ed.A

Copyright© ALPHADENT NV 1/2008

BTW BE 0414.789.321 RPR Kortrijk

ALPHADENT NV

Textielstraat 24

8790 Waregem, Belgium

tel +32 56 629 100

fax +32 56 629 101

e-mail: info@ceka-preciline.com